
kura*ray*





CLEARFIL MAJESTY ES-2

VITA Approved Shade Concept

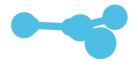

SIMPLY BLENDING INTO NATURE

SIMPLY BLENDING INTO NATURE

CLEARFIL MAJESTY™ ES-2 is a light-cure, radiopaque, universal nano hybrid composite resin, providing accurate colour matching, unsurpassed polishability and excellent physical properties. Unlike traditional composites, CLEARFIL MAJESTY™ ES-2 relies on the initial selection of shade value instead of hue/chroma. This advancement has lead to CLEARFIL MAJESTY™ ES-2 being the only composite to offer VITA Approved Shades and Shade Concept.

Innovative Resin Technology

CLEARFIL MAJESTY™ ES-2 has five Biomimetic characteristics to perfectly simulate natural tooth.



Fluorescence



Light Diffusion



Opacity



Transparency



Opalescence

Versatile Handling

Versatile usage ___

CLEARFIL MAJESTY™ ES-2 is a nano hybrid, restorative composite resin, ideal for anterior and posterior applications. CLEARFIL MAJESTY™ ES-2 is available in a single layer (Classic) and our dual layer (Premium) range.

Enamel paste viscosity

The enamel paste viscosity is optimized for precise placement and sculpting. Adaptation of CLEARFIL MAJESTY™ ES-2 enamel paste to prepared bevels and/or cervical dentine is seamless.

Dentine paste handling _


Enhanced monomer and resin composition delivers outstanding resin integration to the cavity.

Extended working time ____

CLEARFIL MAJESTY™ ES-2 has an extended working time of 4.5 minutes under ambient light enabling ample time for optimal creativity and superior clinical results.

Intuitive Shade Guide

VITA Approved Shade Concept Begin the shade selection process by choosing from 6 value shade tabs for a simplified shade procedure.

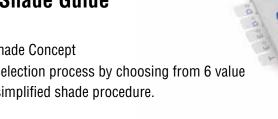

VITA Approved Shades & Shade Concept


Classic: Single shade technique with VITA Approved Shades



Premium: Dual shade technique with our VITA Approved Shade Concept

Biomimetic layering system

Natural dentine replaced with dentine shade and natural enamel replaced with enamel shade

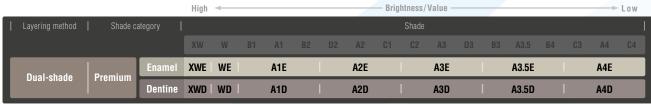

Additional shades: Translucent shades can be used for detailed characterization of the restoration



■ CLEARFIL MAJESTY™ ES-2 Clinical Case

Classic: Single shade technique



■ Classic: Single shade technique with VITA Approved Shades






■ CLEARFIL MAJESTY™ ES-2 Clinical Case

Premium: Dual shade technique



■ Premium: Dual shade technique with our VITA Approved Shade Concept

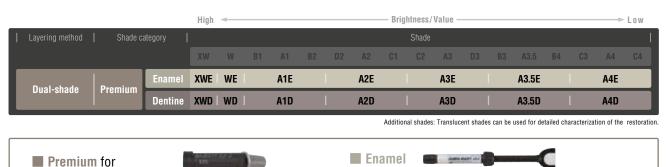



VITA Approved Shade Concept

Select the closest brightness/value tab from 6 tab options on the CLEARFIL MAJESTY™ ES-2 shade guide.

Restoration


- Select the corresponding pastes from the Premium range (Enamel & Dentine). For example, if A1 is selected, A1E and A1D will be used.
- Replace the dentine space with the respective dentine shade and the enamal space with the matching enamel shade.



Expansive Shade Coverage


One CLEARFIL MAJESTY™ ES-2 shade combination covers three VITA Classical shades.

Seeing is believing. The new CLEARFIL MAJESTYTM ES-2 has expansive shade coverage with a condensed shade range. This is possible due to the creation of brightness (value) based shade groups.

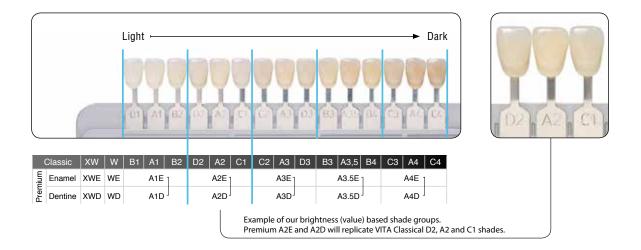

Research: Shade Variance Study

In this study, the centre of three VITA Classical Shade Guide tabs (D2, A2, C1) were restored according to the shade instructions of three manufacturers. With the A2E and A2D shade combination, Kuraray Noritake Dental's CLEARFIL MAJESTY™ ES-2 proved to match not only the A2 shade but also the D2 and C1 shades. All three tabs, restored with CLEARFIL MAJESTY™ ES-2, showed no visible difference to the original VITA tab.

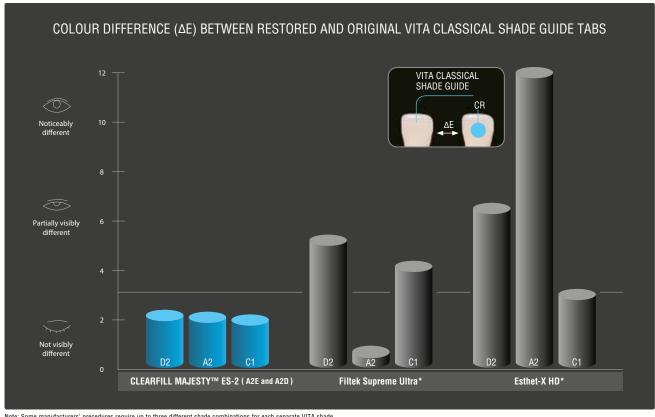


Note: Some manufacturers' precedures require up to three different shade combinations for each separate VITA shade.

CLEARFIL MAJESTY™ ES-2 matches all three VITA shades with only one Premium shade combination (in this case A2E and A2D).


**Not A treating of KIMBADA COLLTD.

* Not a trademark of KURARAY CO., LTD. Source: Dr. Shigemi Nagai, Harvard School of Dental Medicine, United States

■ CLEARFIL MAJESTYTM ES-2 Classic Shades Available in Syringe and PLT refills

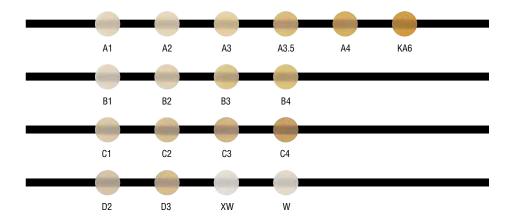

■ CLEARFIL MAJESTYTM ES-2 Premium Shades Available in Syringe and PLT refills

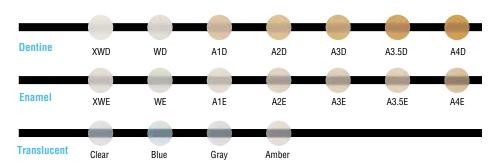

The colours of the actual shade guide may deviate from printed colours.

CLEARFIL MAJESTY™ ES-2

#3010-WD Syringe Professional Kit

Classic A1, A2, A3, A3.5, Premium A1E, A2E, A3E, A3.5E, A1D, A2D, A3D, A3.5D: 1 syringe per shade (3.6 g / 2.0 ml) + Clear, Amber, Blue, Gray 1 syringe per shade (1.8 g / 0.9 ml) + Shade Guide

#3100-WD PLT Professional Kit

Classic A1, A2, A3, A3.5, Premium A1E, A2E, A3E, A3.5E, A1D, A2D, A3D, A3.5D, Clear, Amber, Blue, Gray: 5 tips per shade (0.25 g) + Shade Guide

#3012-WD Syringe Introductory Kit

Classic A2, Premium A2D, A2E: 1 syringe per shade (1.2 g / 0.7 ml)

#3102-WD PLT Introductory Kit

Classic A2, Premium A2D, A2E: 5 tips per shade (0.25 g)

Shade Guide

#3000-WD Shade Guide

All shades

#3001-WD Shade Guide Compact

A1, A2, A3, A3.5, A4, B1, B2, B3, B4, C1, C2, C3, C4, D2, D3, XW, W

PLT Refill Classic

#3110-WD	A1	(20 tips x 0.25 g)	#3140-WD	A1E	(10 tips x 0.25 g)
#3111-WD	A2	(20 tips x 0.25 g)	#3141-WD	A2E	(10 tips x 0.25 g)
#3112-WD	A3	(20 tips x 0.25 g)	#3142-WD	A3E	(10 tips x 0.25 g)
#3113-WD	A3.5	(20 tips x 0.25 g)	#3143-WD	A3.5E	(10 tips x 0.25 g)
#3114-WD	A4	(20 tips x 0.25 g)	#3144-WD	A4E	(10 tips x 0.25 g)
#3115-WD	KA6	(20 tips x 0.25 g)	#3145-WD	XWE	(10 tips x 0.25 g)
#3116-WD	B1	(20 tips x 0.25 g)	#3146-WD	WE	(10 tips x 0.25 g)
#3117-WD	B2	(20 tips x 0.25 g)	#3147-WD	A1D	(10 tips x 0.25 g)
#3118-WD	В3	(20 tips x 0.25 g)	#3148-WD	A2D	(10 tips x 0.25 g)
#3119-WD	B4	(20 tips x 0.25 g)	#3149-WD	A3D	(10 tips x 0.25 g)
#3120-WD	C1	(20 tips x 0.25 g)	#3150-WD	A3.5D	(10 tips x 0.25 g)
#3121-WD	C2	(20 tips x 0.25 g)	#3151-WD	A4D	(10 tips x 0.25 g)
#3122-WD	C3	(20 tips x 0.25 g)	#3152-WD	XWD	(10 tips x 0.25 g)
#3123-WD	C4	(20 tips x 0.25 g)	#3153-WD	WD	(10 tips x 0.25 g)
#3124-WD	D2	(20 tips x 0.25 g)	#3154-WD	Clear	(10 tips x 0.25 g)
#3125-WD	D3	(20 tips x 0.25 g)	#3155-WD	Amber	(10 tips x 0.25 g)
#3126-WD	XW	(20 tips x 0.25 g)	#3156-WD	Blue	(10 tips x 0.25 g)
#3127-WD	W	(20 tins x 0.25 a.)	#3157-WD	Grav	(10 tins x 0.25 a.)

Premium

Syringe Refill

Oyimge	, ,,,	,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,				
Classic			Premium			
#3020-WD	A1	(3.6 g / 2.0 ml)	#3050-WD	A1E	(3.6 g / 2.0 ml)	
#3021-WD	A2	(3.6 g / 2.0 ml)	#3051-WD	A2E	(3.6 g / 2.0 ml)	
#3022-WD	A3	(3.6 g / 2.0 ml)	#3052-WD	A3E	(3.6 g / 2.0 ml)	
#3023-WD	A3.5	(3.6 g / 2.0 ml)	#3053-WD	A3.5E	(3.6 g / 2.0 ml)	
#3024-WD	A4	(3.6 g / 2.0 ml)	#3054-WD	A4E	(3.6 g / 2.0 ml)	
#3025-WD	KA6	(3.6 g / 2.0 ml)	#3055-WD	XWE	(3.6 g / 2.0 ml)	
#3026-WD	B1	(3.6 g / 2.0 ml)	#3056-WD	WE	(3.6 g / 2.0 ml)	
#3027-WD	B2	(3.6 g / 2.0 ml)	#3057-WD	A1D	(3.6 g / 2.0 ml)	
#3028-WD	B3	(3.6 g / 2.0 ml)	#3058-WD	A2D	(3.6 g / 2.0 ml)	
#3029-WD	B4	(3.6 g / 2.0 ml)	#3059-WD	A3D	(3.6 g / 2.0 ml)	
#3030-WD	C1	(3.6 g / 2.0 ml)	#3060-WD	A3.5D	(3.6 g / 2.0 ml)	
#3031-WD	C2	(3.6 g / 2.0 ml)	#3061-WD	A4D	(3.6 g / 2.0 ml)	
#3032-WD	C3	(3.6 g / 2.0 ml)	#3062-WD	XWD	(3.6 g / 2.0 ml)	
#3033-WD	C4	(3.6 g / 2.0 ml)	#3063-WD	WD	(3.6 g / 2.0 ml)	
#3034-WD	D2	(3.6 g / 2.0 ml)	#3064-WD	Clear	(1.8 g / 0.9 ml)	
#3035-WD	D3	(3.6 g / 2.0 ml)	#3065-WD	Amber	(1.8 g / 0.9 ml)	
#3036-WD	XW	(3.6 g / 2.0 ml)	#3066-WD	Blue	(1.8 g / 0.9 ml)	
#3037-WD	W	(3.6 g / 2.0 ml)	#3067-WD	Gray	(1.8 g / 0.9 ml)	

Kuraray Noritake Dental Inc.

In 1926, Kuraray was established in Kurashiki, Japan, where the industrial production of viscose rayon began. Though intensive research and development the company succeeded in building up extensive expertise and technological strength in the fields of polymer chemistry, synthetic chemistry and chemical engineering.

Kuraray entered the field of dental materials in 1973.

The mission of Kuraray was to respond precisely and thoroughly to the needs of the dental profession with reliable products of superior quality. In 1978, Kuraray introduced the first bonding system to the market:

CLEARFIL™ BOND SYSTEM F, the start of the age of adhesive dentistry. At the same time, the company developed the "total-etch technique" for enamel and dentine.

Over four decades, Kuraray Noritake Dental Inc. continues to bring forth innovative, quality products to meet the demands of a continuously changing dental profession.

Ote Center Bldg., 1-1-3, Otemachi, Chiyoda-ku, Tokyo 100-0004 http://www.kuraraynoritake.com/ANZ

[•] Read the Instructions for Use supplied with the products before use.

[•] The specifications and appearance of the product are subject to change without notice.