



HYDRIM[®] G4

Washer-Disinfector






LEADERSHIP

in Infection Control

Instrument washing is the first step in the sterilisation process, which minimizes the physical handling of instruments and the many steps involved with manual cleaning such as pre-soaking, scrubbing and drying instruments. In addition to meeting the requirements of various guidelines and controls, automated washing and disinfection provides a safe, effective and efficient way to clean instruments, allowing staff to focus on the most important thing in the dental office, the patient.

Through a sustained legacy of product innovation and refinement, SciCan developed the HYDR/M[®] Washer-Disinfectors, which provide a suitable solution for nearly every office, whether as a stand-alone or counter-top unit.



Originally founded in 1957, SciCan develops and manufactures highly innovative infection control products that address the need for speed, flexibility, and efficiency in today's environments.

In 1989, SciCan developed and produced an internationally recognized infection control innovation, the STAT/M® Cassette Autoclave. Over the next quarter century, SciCan continued to be at the forefront of technology with the release of the HYDRIM® Instrument Washer-Disinfectors, BRAVO™ Chamber Autoclaves, and G4 Networking Technology.

SciCan is a global leader in infection control specializing in dentistry and ophthalmology in more than 100 countries around the world. The portfolio of products continues to grow to address market needs; innovative products such as the SALUS® Instrument Reprocessing container, high-quality turbines, handpieces and contra angles with the appropriate cleaning and maintenance unit, the STATMATIC™ Smart, and the OPTIM® Cleaner and Disinfectant complete the product range.

In 2007, SciCan opened the European headquarters, located in Germany. Our team consists of highly motivated and technically informed individuals whose network of authorities, infection control specialists and other manufacturers allow us to provide legally compliant and evidence-based information. We work closely with specialists in the dental and medical field, enabling us to provide practical guidance and help drive product improvements that will benefit our customers.

It is very important to me and my team, especially in the time of online sales, to offer you not only the best products but also a service that makes our customers feel well advised and cared for. Please don't hesitate to contact our office for any infection control related questions.

Good luck. I wish you only good days and sufficient time for the essential things.

Stefan Helsing (Managing Director)

HYDRIM G4

Instrument Washer-Disinfector

HYDRIM G4 automated instrument washer-disinfector is the first step in the sterilisation process.



Reliable

Designed in accordance with EN ISO 15883. The washer- disinfector has validated cycles to ensure repeatable performance. HYDRIM meets all EU Standards and RKI Requirements related to washer-disinfectors.



Safe

Automated washer-disinfectors reduce physical handling of instruments, decreasing the risk of exposure and injuries. HYDRIM G4 ensures that every set of instruments receive the same validated cleaning process by dosing the correct amount of detergent with each cycle.



Dry

Active air drying technology with hepa filter quickly dries the instruments and allows them to immediately be transferred to the next step of instrument reprocessing, sterilisation.



Capacity

Various options for baskets and adaptors ensure you have a unit that meets the needs of any practice. HYDRIM C61WD has a 6 full-size cassette capacity, allowing offices to turn over 120 instruments in just 60 minutes – start to dry.*



Compact Design

The HYDRIM C61WD, as a table top unit, is perfect for smaller offices. It can be installed under the counter or as a freestanding unit. The cleaning solution is stored in the unit, so there is no extra storage space required.



G4 Technology

Allows for reliable load release; connects and sends cycle data to your smart devices or computer, providing easily accessible cycle data, maintenance information and service logs. Should you wish, G4 technology can connect your HYDRIM to specialized technicians, which allows for remote support.

*Actual cycle times will vary with incoming water temperature and pressure.

HYDRIM C61WD Specifications



Size	Dimensions: (w x h x d)	598 x 520 x 526 mm
	Depth with door open	829 mm
	Weight	44 kg
Cycle Times*	Rinse and Hold	5 mins
	Heavy Duty Wash with Disinfection + Dry	42 mins. + 10 mins.
Water Consumption Per Cycle†	Rinse and Hold	11 litres
	Heavy Duty Wash & Disinfection	30 litres
Electrical Connection‡	Electrical Rating	208-240V or 200-230V, 50 Hz or 60 Hz, 12A, single phase
Water and Drain Connections	Hot and Cold Water Tap Connection	3/4" Garden hose, male
	Recommended Hot Water Temperature	60°C
	Cold Water Temperature	< 30°C
	Inlet Water Pressure	1 - 10 bar
	Drain	2 cm

* Actual cycle times will vary with incoming water temperature and pressure.

† Water consumption includes drying

‡ Electrical requirements are model specific

HYDRIM C61WD Rack Configurations

Configuration 1



3 Cassette /
2 Basket Rack
01-113253



Handpiece
system LCS
01-113625



Vertical
Instrument Rack
01-113255



Bar for Scaler Tip
with Torque Tool
01-114841

Additional Accessories

Adapter for
Handpieces
01-114845

Metal Plate
Holder for Bars
01-114846

Configuration 2




Handpiece
system LCS
01-113625



Vertical
Instrument Rack
01-113255



Bar for Scaler Tip
with Torque Tool
01-114841

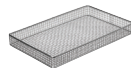

STATIM 2000
Basket
01-107240



Basket Rack
01-113257



Lower Trolley
Basket
01-114942

Additional Accessories

Adapter for
Handpieces
01-114845

Metal Plate
Holder for Bars
01-114846

Configuration 3



3 Cassette /
2 Basket Rack
01-113253



Handpiece
system LCS
01-113625



Bar for Scaler Tip
with Torque Tool
01-114841



Multi-Use Bar
Accessory
01-114842

Additional Accessories

Adapter for
Handpieces
01-114845


Metal Plate
Holder for Bars
01-114846

Adapter for Surgical
Suction Tubes
01-114844


Configuration 4



Full Basket
01-113545

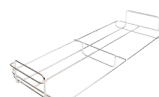

Basket Rack
01-113257

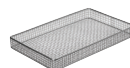

STATIM 2000
Basket
01-107240



STATIM 5000
Basket
01-107241

More loading samples available upon request.

Cleaning

Cleaning is one of the most important steps in the instrument reprocessing procedure. In order to achieve effective disinfection and sterilisation, the instruments need to properly be cleaned first. Washer-Disinfectors have been proven to be the most safe and effective way to clean and disinfect instruments, which is supported by many regulating authorities and dentists worldwide.



Made of robust stainless steel, with a comprehensive range of accessories, and solutions for almost every office. No extra detergent space is required as the patented cleaning solution is stored inside the unit. The HYDRIM washer-disinfector is user-friendly, safe and effective, allowing for thorough cleaning as well as thermal disinfection at 93°C.

By default all HYDRIM units offer active air drying with a hepa filter, protecting instruments from recontamination and corrosion.

Conformity

HYDRIM C61 WD G4 are produced in conformity with the following directives and standards:

93/42/EEC	Medical device class IIa (European Directive for Medical Devices)	EN 1717	Protection against Pollution of Potable Water
EN ISO 15883-1/-2	Washer-Disinfectors	EN ISO 13485:2016	Medical devices – Quality management systems – Requirements for regulatory purposes
EN 61010-1 /-2-040	Safety Requirements for Electrical Equipment for Measurement, Control and Laboratory Use	EN ISO 14971:2012	Application of risk management to medical devices
EN 61326-1	Electrical equipment for measurement, control and laboratory use – EMC requirements	KIWA certified	

G4 Technology

As special feature of HYDRIM units is the innovative G4 technology. A large touch screen provides in-depth cycle info and allows users to view current operation.

With an integrated Ethernet port, HYDRIM can connect to the office network and send cycle data to your smart devices or computer, providing easily accessible cycle data, maintenance information and service logs.

To help optimize office processes, the HYDRIM's can be configured to email cycle information or error messages directly to office staff or service technicians. Cycle data is automatically saved internally, on a USB and can also be emailed or printed with a network printer, a feature that has been specifically designed for meeting the growing legal demand of physical documentation.

The HYDRIM G4 units have a traceability features that requires a user code when releasing a load.



HYDROS™ - Water Treatment System

for use with washer-disinfectors

Bad water quality can lead to discoloration and calcium deposits in the machine. Connecting the washer-disinfector to a R/O water treatment system is therefore recommended to achieve optimal cleaning and support the longevity of your machine.

Demineralized water (R/O water) is also recommended in the RKI guideline from 2012 (Hygiene requirements for the reprocessing of medical devices) for the final rinse in washer-disinfectors.

HYDROS has a special multi-layered filter that removes several unwanted components of the water, step by step.



User Friendly Operation

Easy exchange of the filter cartridge



Space Saving Installation

Horizontal and vertical installations possible



Easy To Use

Suited for units with hot water connection
(up to 60°C)



Technical Specifications

Capacity at a carbonate hardness of 10° dH *	5000 l
Operating pressure	6 bar
Water intake temperature	4 – 60°C
Flow at 1 bar pressure loss	850 l/h
Nominal flow	300l/h
Pressure loss at nominal flow	0.45 bar
Dimensions (width/depth/height)	288/255/550 mm
Weight (dry/wet)	18 / 24 kg
Connections (input/output)	G1" / G3/4"

* The stated capacity is a guide value only, that may vary by +/- 20% depending on local water quality and machine type.

SciCan's Full Spectrum Solutions

Helping you get the most out of your sterilisation area and the equipment you put in it.





Manufactured by:

SciCan Ltd.

1440 Don Mills Road, Toronto,
ON, Canada M3B 3P9

Tel: 416-445-1600

Fax: 416-445-2727

Toll free: 1-800-667-7733

customerservice@scican.com



EU Representative & Distributed by:

SciCan GmbH

Wangener Strasse 78

88299 Leutkirch

Deutschland

Tel.: +49 (0)7561 98343 - 0

Fax: +49 (0)7561 98343 - 699

customerservice_eu@scican.com

Office for Swiss Customers:

SciCan Medtech AG

Alpenstrasse 14

6300 Zug

Switzerland

Tel: +41 (0) 41-727 7027

Fax: +41 (0) 41-727 7029