

BL S-Kondenser

Warm gutta percha obturation techniques usually require condensation with vertical or lateral forces for an effective 3-dimensional fill. Essential to this technique are precision condensers which can condense thermoplasticized gutta percha in complex anatomy, including extremely curved canals and ovoid shaped canals.

The BL S-Kondensers have been designed specifically for warm gutta percha techniques with these important design features:

- A .04 taper NiTi end for apical condensation as well as a larger stainless steel end for middle and coronal third condensation
- A size 35 for smaller canal diameters
- NiTi end navigates even extremely curved canals and is manufactured with the highest grade Nickel Titanium
- Light-weight anodized aluminium handles
- Ergonomically placed finger rests for efficient condensing
- Laser etched markings on the NiTi end to identify the depth of the condenser



Ordering Codes:

- BCON35** BL S-Kondenser 35/70
BCON40 BL S-Kondenser 40/80
BCON50 BL S-Kondenser 50/100
BCON60 BL S-Kondenser 60/120
BCONSET BL S-Kondenser Set (35, 40, 50, 60)

SuperEndo™ α^2 & β

Cordless Obturation Device

Compact, Convenient, Smart and Versatile,
 Consistent heat, 3D-Obturation



SuperEndo™ α²

Cordless Obturation Device



SuperEndo™ β

Cordless Obturation Device



Introducing the Alpha II to enhance the success of conventional obturation techniques. Simple to use, this innovative device is compact, convenient and versatile. Alpha II's unique design offers 3-dimensional obturation, temperature control, as well as the convenience and ergonomics expected of a cordless product.


- Lightest device on the market
- Precision control of heat
- A variety of pluggers and instruments to suit any tooth and warm GP obturation technique
- Variable positions for easy access into any canal
- Very long battery life
- Quickly reaches desired temperature setting



Plugger Stand
 Storage for B&L BL Pluggers

Ordering Codes:

ASK	Alpha II Starter Kit
BAB	Alpha II Replacement Battery
HTT	Hot Testing Tip
BPSS	Plugger Protector Sleeve
BP3004	30/04 Taper Plugger
BP3504	35/04 Taper Plugger
BP4004	40/04 Taper Plugger
BP4504	45/04 Taper Plugger
BP5506	55/06 Taper Plugger
BP5508	55/08 Taper Plugger
BP5508L	55/08 Taper Plugger Long Type
BP5510	55/10 Taper Plugger
BP6012	60/12 Taper Plugger
SPRDN	Endo Spreader (Narrow)
SPRDS	Endo Spreader (Standard)

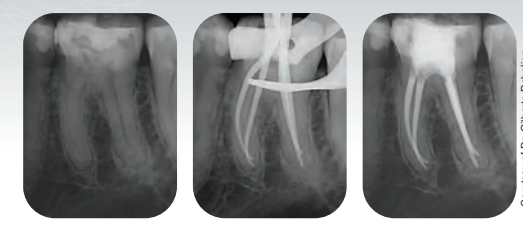

Obtured with warm gutta percha technique using Alpha II & Beta

Courtesy of Dr. Gilberto Debellin

CORDLESS. COMPACT.

With a battery that lasts days on a single charge*, the BETA is the ultimate in portability. The small footprint of the BETA charge cradle makes the unit both convenient and compact.

ERGONOMIC. PRACTICAL.

Lightweight and balanced for comfort, the BETA packs its punch in durability, withstanding the wear and tear of repeated use in the busiest of endodontic practices. This virtually maintenance-free unit also features a well-insulated heating chamber that keeps the gun surface cool to the touch for your safety as well as your patient's.

* On full charge, the BETA delivers 4 hours of continuous operation.

B&L BL Pluggers

- Durable, high-grade stainless steel
- Precise temperature control
- Cools very quickly once heat is turned off
- Largest selection of sizes



DESIGNED FOR BETTER DENTISTRY.

With the BETA in hand, you are in control. At any time during use you may adjust the temperature and/or your access angle to maintain ideal obturation flow.

4 TEMPERATURE SETTINGS – With four accurate temperature pre-sets and an instant on/off at your fingertips, you control the consistency and temperature of the gutta percha or Resilon at all times during obturation, without wasting battery power.



360° SWIVEL NEEDLE – The primary clinical challenge with all injectable obturation procedures is the angle of entry to the canal, especially on maxillary cases. The BETA's patented 360° swivel needle permits the clinician to access the canal at any angle with ease. No more awkward gun or hand positions, no more stopping midway through the procedure to adjust the needle angle, and no more needle breakage.

A PERFECT COMPLEMENT TO ALPHA2 – For the ultimate in warm vertical obturation efficiency, pair ALPHA A2 Cordless Heat Source with BETA Cordless Obturation Gun.



Gutta Percha Pellets
 (available in soft, regular and hard)



Needles
 (available in 21mm, 22mm, 24mm, 28mm & 50mm lengths and 20, 23 and 25 gauges)

Ordering Codes:

BSK	Beta Starter Kit
BTP	Beta Thermal Protectors
BPA	Beta Plugger Assembly
BNPS	Beta Needle Protector Sleeve
BGP####	GP Pellet (SOFT / REG / HARD)
BSN2123	21mm Needle, 23Gauge
BSN2220	22mm Needle, 20Gauge
BSN2420	24mm Needle, 20Gauge
BSN2423	24mm Needle, 23Gauge
BSN2425	24mm Needle, 25Gauge
BSN2625	26mm Needle, 25Gauge
BSN2820	28mm Needle, 20Gauge
BSN2823	28mm Needle, 23Gauge
BSN5020	50mm Needle, 20Gauge - VET
BMT	Needle Bender / Wrench