Good care pays off.

The cost-efficient solution for optimum care and long-term maintenance of value.











KaVo. Dental Excellence.

Perfect quadruple care. Fast and comfortable.

Very easy to use.

Normally, there is little available time in a busy practice, to provide adequate care for high-quality handpieces. Little time for you and your practice team. Thus, it is very important to rely on a sophisticated care system, with intuitive, fast and easy handling. Mount handpieces. Close front panel. Start. Four at a time. Finished in no time at all.

Four handpieces in one minute.

Just by the press of a button. Simple and fully automated. KaVo QUATTROcare PLUS, offers maximum speed and user-friendliness. The KaVo QUATTROcare PLUS is characterised by easy handling. Mount up to four handpieces on the four couplings, close the front panel and start the care program. The complete care procedure for four handpieces, only takes one minute. Then the QUATTROcare PLUS, can be immediately used for the next care cycle.



It's as simple as that.

Optimum care quality, for up to four handpieces:

The procedural steps with KaVo QUATTROcare PLUS.

- Mount handpieces. One MULTIflex treatment coupling, as well as three INTRAmatic treatment couplings, are supplied as standard for mounting of handpieces.
- Close front panel and you are ready to go.
- Press start button. The fully automated, internal cleaning and care, of up to four handpieces, only takes approx. 1 minute.
- Finished. Remove the handpieces.











- No operating errors, thanks to easy and safe handling
- Care safety, thanks to automatic program sequence
- Renewed, availability of handpieces, after only one minute, thanks to the short cycle-time

KaVo QUATTROcare PLUS: simple - fast - efficient.





Optimum care.

The difference, as always, is in the details.

The rapid expansion of the foaming, cleaning agent in the KaVo QUATTROcare PLUS "Care Spray", dissolves inorganic material within the handpiece and then compressed-air blows out the residue, thereby providing the basis for optimum handpiece care and maintenance.

Although not essential, previous application of both KaVo CLEAN and KaVo DRYspray is recommended, to definitively remove proteins from within the handpiece, prior to lubrication. The speed of the KaVo QUATTROcare PLUS, not only expedites your routine practice procedures, but also provides perfect handpiece care. This has a very positive effect on the service-life of your handpieces.

The automatic, metered dosing, ensures that each handpiece is subjected to maintenance with the prescribed amount of oil, corresponding to its type. Thus, unnecessary use of care oil is avoided and optimum handpiece performance is assured. KaVo QUATTROcare PLUS is universally compatible. This means that you can use KaVo QUATTROcare PLUS, not only for the care of your high-quality KaVo handpieces, but also for the care of handpieces from other established manufacturers.

EXCELLENT.

Consequent efficiency of care.

Longevity is not only a question of product quality, but also of regular, thorough maintenance. Thanks to the KaVo QUATTROcare maintenance system, KaVo sets standards of efficiency for hand-piece care. The innovative, KaVo QUATTROcare PLUS maintenance system, offers incomparable thoroughness, economy and easy-of-use. Perfect for every handpiece used at your practice, every day.

KaVo QUATTROcare PLUS is equipped with the innovative, purging function, as standard. The handpiece transmission and spray channels, are additionally supplied with compressed-air and purged, so that any residual oil or detritus is removed and all handpieces therefore, receive truly comprehensive care.



Integrated, chuck-care program.

Once-per-week, cleaning and care of handpiece chucks, maintains their gripping force, reliability and long service-life

- 1. Mount chuck-care adapter onto MULTIflex coupling
- 2. Place chuck into care adapter
- 3. Press program key for chuck-care









The KaVo QUATTROcare effect.
Powerful, foaming expansion of the KaVo QUATTROcare
PLUS "Care Spray" within the
handpiece, dissolves and removes
organic residues.



SUPERtorque LUX 650 B - care with KaVo QUATTROcare maintenance system. Seven years, in daily use at the practice - without failure or replacement of the cartridge. Always optimum care - always ready-to-use.

Efficiency at the practice. KaVo

Care and efficiency – main features and advantages for you:

- Optimum care, thanks to KaVo QUATTROcare PLUS "Care Spray"
- High efficiency, thanks to being perfectly adapted to daily practice requirements

KaVo QUATTROcare PLUS: optimum care – highly efficient.



Favourable purchase cost. Economical.



Benefits are obvious.

Do you invest in top-quality handpieces? KaVo QUATTROcare PLUS ensures that your instruments maintain precision and value. Continuous cost-savings and optimum performance, are only possible with perfect care and hygienically Reprocessed handpieces.

Automatic, metered dosing and the automatic, outlet-specific, handpiece detection system of the KaVo QUATTROcare PLUS, ensure perfect, precise and economic consumption of KaVo QUATTROcare PLUS "Care Spray". Unnecessary use and correspondingly increased costs cannot occur.



Cost-efficient care – main features and advantages for you:

- Only one-minute care cycle for four handpieces
- Program-controlled, automated care process
- Optimum lubrication, thanks to automatic, outlet-specific, metered dosing
- Easy handling prevents erroneous care
- Compressed-air blasting of transmission system and spray channels, using purging function
- No leakage of excess, residual oil from handpiece, after processing
- Basic device, complete with four treatment couplings (1 x MULTIflex / 3 x INTRAmatic)
- Universal, handpiece compatibility (adapters for all established brands available)
- Modern and ergonomic design

KaVo QUATTROcare PLUS: simple - versatile - economical.





Mat.-Nr. 1.004.1434 03/11 en We reserve the right to make technical modifications. Slight colour differences are due to the printing process. © Copyright KaVo Dental GmbH.

Technical Specifications

Dimensions:


Width 35 cm
Depth 28 cm
Height 38 cm

Technical data:

Weight	10 kg
Mains voltage	100 – 240 V
Mains frequency	50 – 60 Hz
Power consumption	approx. 24 W
Compressed-air pressure	4 – 6 bar
Compressed-air consumption	approx. 50 NI/min

Installation requirement:

- good room ventilation
- no direct exposure to sunlight

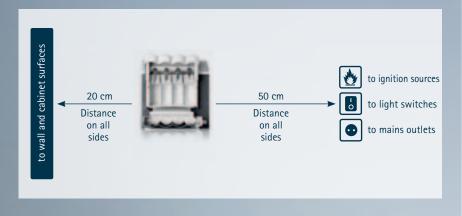

Accessories

Product	Mat. no.
INTRAmatic treatment coupling	0.411.7981
MULTIflex treatment coupling	0.411.7991
INTRA Head treatment coupling	0.411.7941
COMFORTdrive treatment coupling	1.005.1707
SIRONA TE treatment coupling	1.000.1003
SIRONA T1 turbine, treatment coupling	1.000.7156
BIEN-AIR treatment coupling	1.001.7860
W&H treatment coupling	1.001.7888

Product	Mat. no.
YOSHIDA turbine with light, treatment coupling	1.000.7403
YOSHIDA turbine without light, treatment coupling	1.000.7585
MIDWEST Low Speed, treatment coupling	1.000.8459
NSK turbine treatment coupling	1.000.8786
OSADA turbine treatment coupling	1.000.8858
Borden 2-hole treatment coupling	1.002.1217
QUATTROcare PLUS Care Spray (package with x6 can	s) 1.005.4525

