




# Fusion™ Anterior Matrix System

**The Fusion Anterior Matrix System from Garrison Dental Solutions allows you to restore many types of anterior lesions with Confidence!**

**The Social Six never had it so good!**


Fusion Anterior Matrix bands help you recreate ideal interproximal anatomy quickly and confidently:

- Only 0.0015" - significantly thinner than strips to facilitate better interproximal contact.
- Firm stainless steel resists deformation both during placement and while compressing composite into the preparation.
- Far superior to strips when restoring large or deep carious lesions.

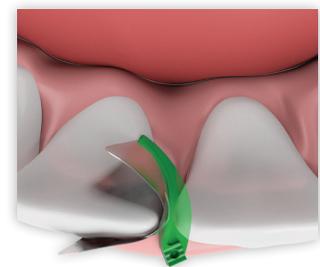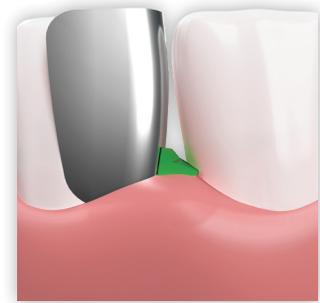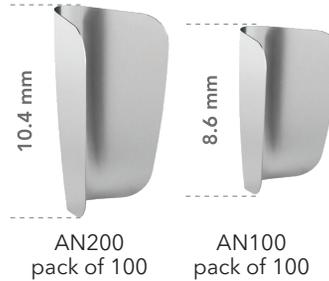

The clever features of the Fusion Anterior Wedge make it a 'must have' for all anterior restorations:

- Radical curvature ensures a firm seal at the cervical margin from facial to lingual.
- Tapered separator strut allows you to adjust the amount of pressure exerted to control the interproximal contact.
- The cervical interface is half above and half below the sulcus for maximum adaptation with deeper restorations.

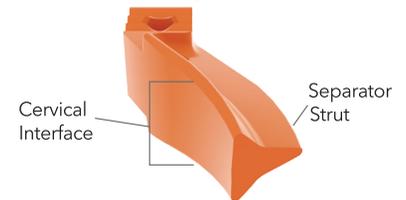

**ANK01 Contains:**  
 50 short bands  
 50 tall bands  
 25 extra small wedges  
 25 small wedges  
 25 medium wedges  
 25 large wedges

Try the all-new **Fusion™ Anterior Matrix System** today!



Refills	Item	Qty.	Price
X-Small Yellow	ANYL	100	
Small Blue	ANBL	100	
Medium Orange	ANOR	100	
Large Green	ANGR	100	

Products available through these authorized dealers.

AUSTRALIA



NEW ZEALAND



www.henryschein.co.nz  
 Orders: 0800 808 855

THE LEADER IN MATRIX SYSTEMS

© 2021 Garrison Dental Solutions, LLC



Everything you order from Garrison Dental Solutions is backed by our risk-free 60-day money-back.



ADAU421 ID