




Real Easy. Real Tight. Real Fast!

For all posterior restorations!

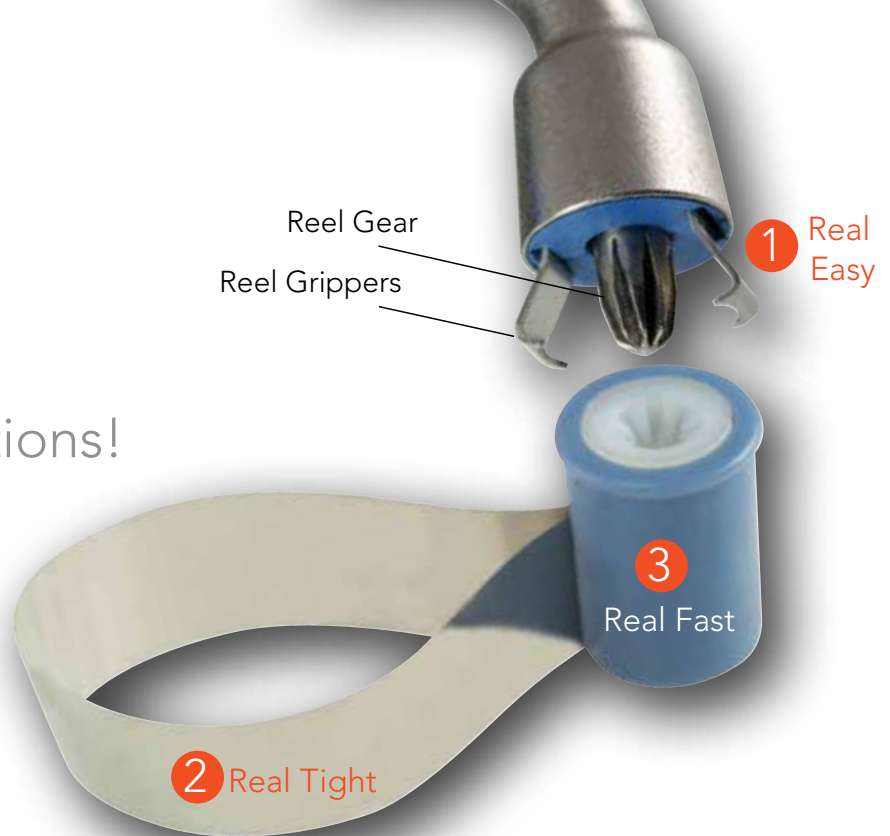

Real Easy!

The reels pop into the handle without any fuss! Press the reel into the handle and retract the grippers to lock the reel in place. Plus, if you use the regular matrices (yellow hubs) they can be easily pushed through contacts.

Real Tight!

Tight contacts and excellent contour are what Reel Matrix is all about. These thin matrices are pre-contoured in three dimensions and when combined with anatomical wedges and separator rings, the result is amazingly tight contacts virtually flash-free at the natural height of contour.

Real Fast!

Pre-loaded reels mean you can just lock them in and go. The retainerless design greatly speeds up quadrant work while improving visibility and patient comfort.

Margin Elevation Band restores deep caries without the trauma.

Margin elevation requires a matrix that will properly adapt and isolate the subgingival area. This is not possible with traditional matrix materials (Figure A). The narrow Margin Elevation Band, with its enhanced conical design, is able to slide down subgingivally and properly seal this area (Figure B).



NEW!
Included in kit!

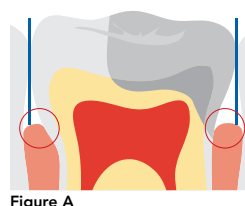


Figure A

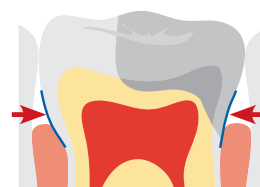


Figure B



Garrison
Dental Solutions

ReelMatrix

Real Easy. Real Tight. Real Fast!

NEW!



RM500H
Margin elevation



RM100D
5mm deadsoft



RM100H
5mm regular



RM100B
5mm blue



RM200D
6.3mm deadsoft



RM200H
6.3mm regular



RM200B
6.3mm blue

Margin elevation band is .002" thin.

Super-thin dead-soft bands (white hubs) are .0012" thin.

Regular bands (yellow hubs) are .0015" thin and can be placed through existing contacts.

A Blue View cure-through option at .0029" thin is available if you prefer trans-enamel curing.



RMR200
6.7mm



RMR100
5.9mm

If there is a band style you really love, reel it up yourself! ReelMatrix fits most standard Tofflemire matrices, including Garrison Slick Bands™. We recommend you try them!



Kit Description

Kit
RMK05
Contains 80 matrices
(20 RM500H, 10 ea. RM100D,
RM200D, RM100H, RM200H,
RM100B, RM200B)
1 tensioning instrument
1 PerForm™ instruments (set of 2)



Refills	Description	Qty.
RMH-100	Tensioning instrument	
RM500H	Margin elevation band	50
RM100D	5mm dead-soft matrices	50
RM200D	6.3mm dead-soft matrices	50
RM100H	5mm regular matrices	50
RM200H	6.3mm regular matrices	50
RM100B	5mm blue matrices	50
RM200B	6.3mm blue matrices	50
RMR100	5.9mm empty reel	100
RMR200	6.7mm empty reel	100



For more information go to www.garrisdental.com. Be sure to check out the website to see the full line of innovative Garrison products and our 6-month money-back guarantee.

Garrison
Dental Solutions

150 DeWitt Lane
Spring Lake, MI USA
49456

To order:
Call toll free **888.437.0032**
or online at www.garrisdental.com



Mention Code
IN0216RM