

Your practice is *our* inspiration.™

NX3 | Perfection in permanent cements



NX3

Defining a new generation in permanent cements



NX3 is a permanent resin cement system with an innovative chemistry for unmatched esthetics, adhesion and great versatility. Delivery systems include an automix syringe for dual-cure indications and a light-cure cement for unlimited work time. NX3 cement can be used for all indirect applications, including veneers.


Inspired by:

Color stability, ease-of-use and cleanup, color match and optimum retention are some of the attributes necessary when choosing a cement—NX3 met all of these expectations.

— Mitch A Conditt, DDS

Research

Most self-cure or dual-cure resin cements use a benzoyl peroxide and tertiary amine pair to initiate polymerization and curing, resulting in discoloration over time due to incompatibility with self-etch adhesives. Now with an amine-free redox initiator system and optimized resin matrix, NX3 is the first truly color-stable adhesive resin cement with enhanced bond compatibility for both total-etch and self-etch adhesives.

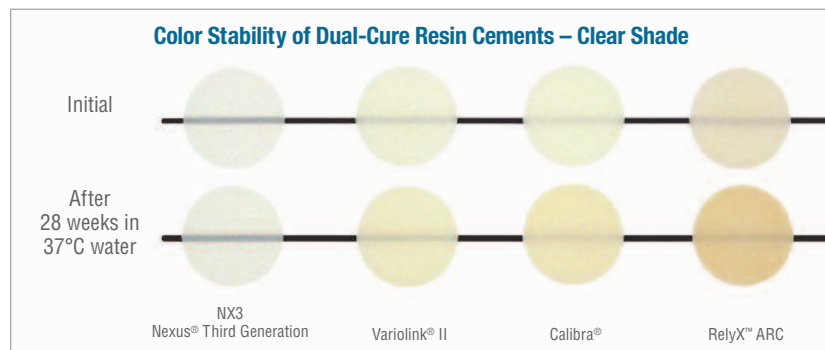


Trademarks are property of their respective owners.

Features

Superior bonding. Bonds to all substrates. Excellent adhesion to dentin, enamel, CAD/CAM blocks, ceramic, porcelain, resin and metal.

Self-etch or total-etch. Bonding protocol compatibility. No dual-cure activator required.

True color stability. Long-term esthetics for both dual-cure and light-cure cements.

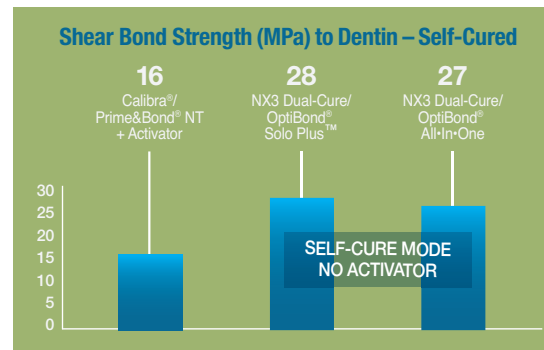

Bond strength study conducted by Dr. James Dunn of Loma Linda University.

Ordering

33642B NX3 Nexus® Third Generation Intro Kit w/risk-free trial offer

Contains 1 Syringe Dual-Cure (5 g) of each shade: clear, white and yellow; 1 Syringe Light-Cure (1.8 g) of each shade: clear, white and yellow; 1 Syringe Try-In Gel (3 g) of each shade: clear, white and yellow; 1 Bottle Silane Primer (5 ml); 30 Automix Tips; Technique Guide; Directions for Use

33682 NX3 Light-Cure Kit

Contains 1 Syringe Light-Cure (1.8 g) of each shade: clear and white; 1 Syringe Try-In Gel (3 g) of each shade: clear and white; 10 OptiBond Solo Plus Unidose Pouches; 1 Bottle Silane Primer (5 ml); 10 Applicator Tips; 1 Syringe Gel Etchant (3 g); 10 Gel Etchant Syringe Tips; Directions for Use

NX3 Automix Dual-Cure Syringe Refills

- 33643 Clear
- 33644 White
- 33645 Yellow
- 33646 Bleach
- 33647 White Opaque

NX3 Light-Cure Syringe Refills

- 33648 Clear
- 33649 White
- 33650 Yellow
- 33651 Bleach
- 33652 White Opaque

NX3 Try-In Gel Syringe Refills

- 33656 Clear
- 33657 White
- 33658 Yellow
- 33659 White Opaque
- 33660 Bleach
- 33655 50 Automix Tips