



TRANSLUCENT



FLEXIBLE




PRECISE


### The advantage is clearly visible.

- **Translucent disc.** Enhances work surface visibility.
- **Flush-mounted mandrel.** Protects neighboring teeth and tissue.
- **Identifiable color-coding system.** Color indicates grit type and abrasive side.
- **Three-step polishing system.** Achieves natural high gloss with fewer steps.<sup>1</sup>
- **Contouring disc.** Extra-Coarse grit available for contouring prior to finishing and polishing.

PRECISE

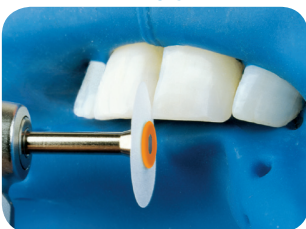


FLEXIBLE

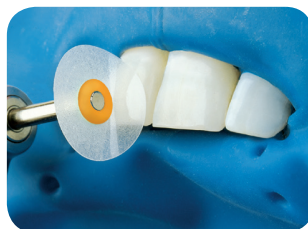

TRANSPARENT



SAFE

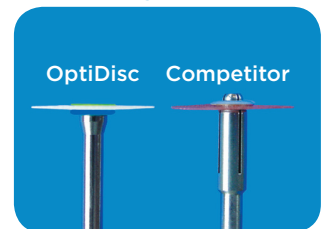


1. Compared to 5 step polishing systems such as 3M Sof-Lex

## OptiDisc Finishing & Polishing Discs

The OptiDisc finishing and polishing system gives your restorations a final polish similar to that of natural dentition. With a three-grit system instead of the traditional four, you can achieve a higher surface gloss with fewer steps than with competitive products. OptiDiscs are translucent for enhanced working surface visibility. Kerr also offers OptiDisc Extra-Coarse for contouring composite restorative materials prior to finishing and polishing.

The color-coded discs clearly identify the grit size. Three disc diameters are available: 9.6 mm, 12.6 mm and 15.9 mm.


## 3-Step OptiDisc™ System vs. 4-Step 3M Sof-Lex™ XT

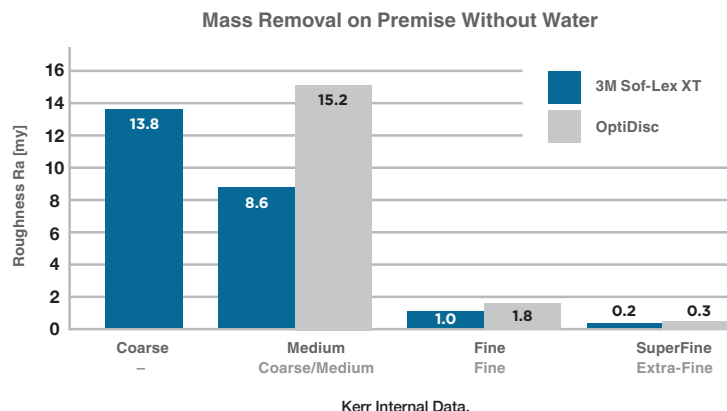

OptiDisc achieves a higher surface polish in just three steps compared with other four-step systems. The charts above demonstrate that the Coarse/Medium OptiDisc has higher mass removal properties than the Coarse Sof-Lex XT, yet inflicts less surface roughness on a restoration, resulting in an easier-to-achieve, smoother final polish.


### Gross Reduction

#### Extra-Coarse

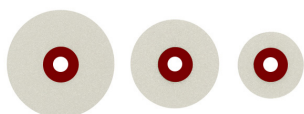


### Polishing

#### Coarse/Medium

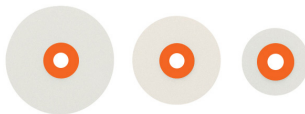

#### Fine

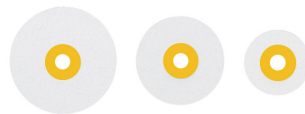

#### Extra-Fine



## Other Features.

- **Long-lasting.** Made of durable polyester, impregnated with aluminum oxide particles.
- **Precise.** Flush-mounted mandrel and patented retention system ensure precise operation, protecting neighboring teeth and tissue.
- **No mess.** Pumice or slurry not required. Use with or without water.

## Ordering

### 33475 OptiDisc Assorted Kit

Contains 30 discs of each:  
 9.6 mm and 12.6 mm – Extra-Fine  
 9.6 mm and 12.6 mm – Fine  
 9.6 mm and 12.6 mm – Coarse/Medium  
 9.6 mm and 12.6 mm – Extra-Coarse 3 Mandrels

### OptiDisc Refills

Contains 30 discs of each:	33482 12.6 mm Coarse/Medium
33476 9.6 mm Extra-Fine	33483 12.6 mm Extra-Coarse
33477 9.6 mm Fine	33484 15.9 mm Extra-Fine
33478 9.6 mm Coarse/Medium	33485 15.9 mm Fine
33479 9.6 mm Extra-Coarse	33486 15.9 mm Coarse/Medium
33480 12.6 mm Extra-Fine	33487 15.9 mm Extra-Coarse
33481 12.6 mm Fine	

