

**THE POWER, SPEED,
AND EFFICIENCY
YOU'VE BEEN WAITING FOR.**



Maxima Advanced MA210
ULTRASONIC CLEANER
BACKED BY A 3-YEAR WARRANTY!

 **HENRY SCHEIN®**
BRAND PRODUCTS

Engineered for powerful cleaning in as little as five minutes*

Now introducing Henry Schein's newest Ultrasonic Cleaner: Maxima Advanced MA210. Our quietest unit to date! With uncompromising performance and design, our compact unit has the **POWER**, **SPEED**, and **EFFICIENCY** you'll want for your practice.



BACKED BY THREE-YEAR WARRANTY



**INDUSTRY
LEADING
CLEAN TIME**



**POWERFUL
CLEANING
TECHNOLOGY**



**QUIET
OPERATION
(69dB)**



**COMPACT SIZE FITS
COMFORTABLY IN
TIGHT SPACES**



**BUILT TO LAST
WITH DURABLE
ALUMINUM**

*Before starting with routine cleaning, test runs should be conducted to establish proper cleaning time. With a standard batch of instruments, run cleaning at 5 minutes then check the condition of the instruments. Add more time until satisfactory result is achieved. Repeat the test a couple of more times to establish a reliable, proper and standard time. Typical cleaning time is 5 minutes for one layer of loose instruments or cassette.



FRONT VIEW



BACK VIEW

DEGAS DIAL AND TIMER

1

to prepare for maximum ultrasonic cleaning performance (60, 120 or 180 sec)

EASY-TO-READ

2

display with bright numbers and icons

START & STOP DIAL/SET TIME

3

for reliable and accurate timer settings (5, 7:30, 10, 12:30, 15, 20 mins.)

QUICK ACCESS

4

to a large 8 Litre capacity tank for easy loading and removal of instruments

SMOOTH CONTOURED EDGES

5

designed with a sleek style in mind

BUILT-TO-LAST

6

durable aluminum panelling

COMPACT SIZE

7

22.9cm height fits comfortably under wall cabinets in your reprocessing area

EASY DRAINAGE ACCESS

8

located on back of the unit so it can be placed where it works best in your reprocessing area

HANDY LID STORAGE

9

located on rear panel to keep your workspace open

ALUMINUM HANDLES

10

for easy grip while transporting machine

Maxima Advanced MA210

Ultrasonic Cleaner HS-570-3794

Large 8 litre tank fits
2 large cassettes or up
to 4 small cassettes

Three stack transducer
technology for power-
ful cleaning



INCLUDED IN BOX: Ultrasonic Cleaning Unit | Stainless Steel Basket Drainage Tube (99cm Long) | Power Cord | Instructions for Use | Quickstart Guide


Features

Tank Capacity	8 Litres
Tank Size	33cm x 23.9cm x 11.7cm
Longest Item Fits Inside Tank	36.6cm
Timer Settings	5 min 7.5 min 10 min 12.5 min 15 min 20 min

Degas Timer Settings	60 sec 120 sec 180 sec
Ultrasonic Frequency	48 KHz $\pm$ 3KHz
Net Weight	14.3kg
Unit Size	46.5cm x 35cm x 22.9cm

Ready to order?

Product Code	Product Name
HS-570-3794	Maxima Advanced Ultrasonic Cleaner MA210
Spare Parts	
HS-570-3824	Power Cord
HS-570-3805	Drain Hose with Elbow
HS-570-3806	Stainless Steel Basket
HS-570-3803	Large Cleaning Chamber
HS-570-3804	Small Cleaning Chamber
HS-570-3802	300mL Beakers Set of 2

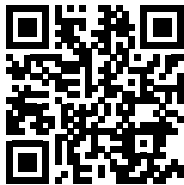

 **HENRY SCHEIN®**
BRAND PRODUCTS

Shop Exclusively at:

HENRYSCHEIN.CO.NZ

243-249 Bush Road, Rosedale, Auckland, New Zealand, 0632
info@henryschein.co.nz | 0800 808 855