


**roeko**

## **Surgitip** ASPIRATOR TIPS

The product family for precise aspiration in surgery and endodontics

**COLTENE**

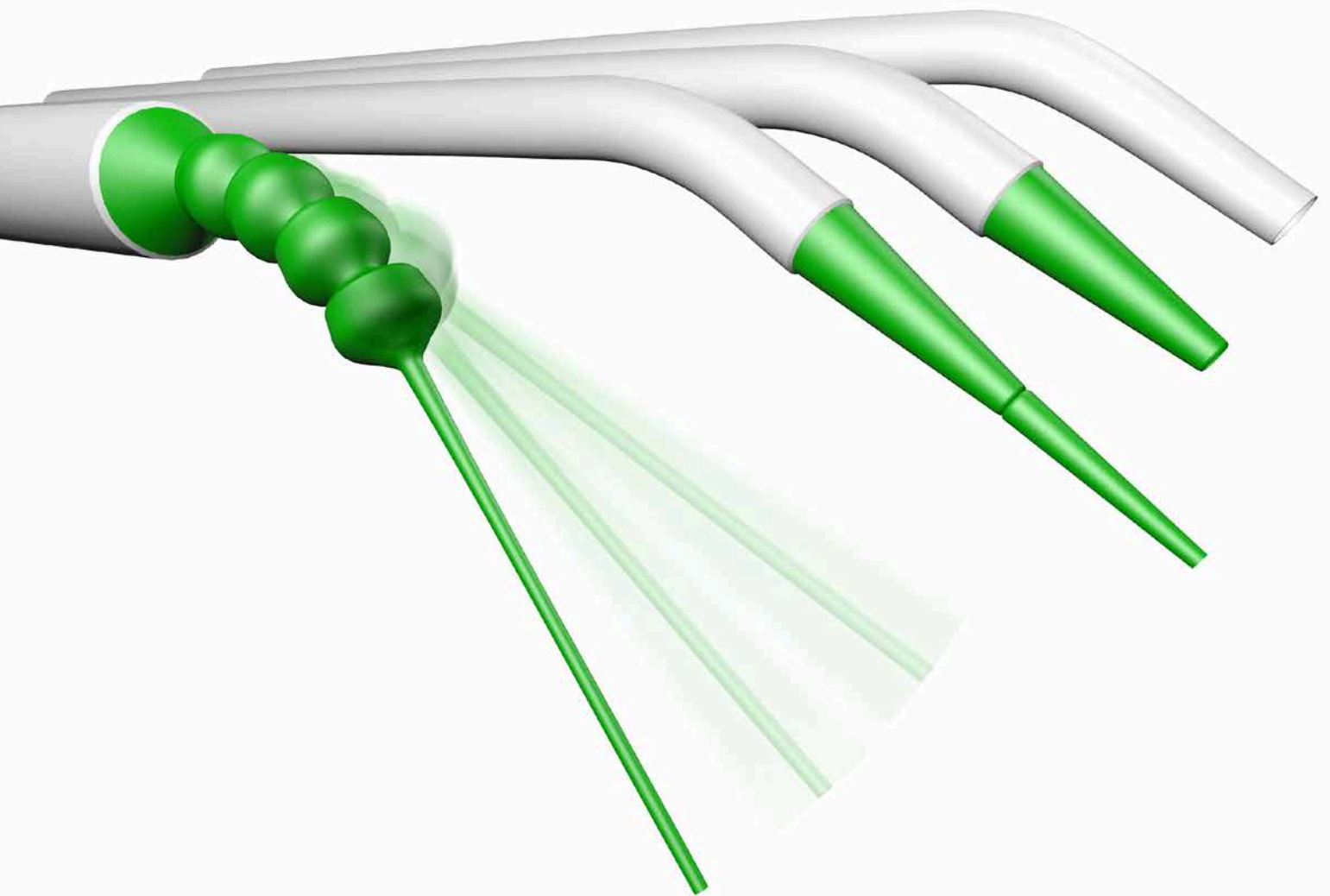

# roeko

## **Surgitip**

Surgitip aspirator tips are proven aids in drying the oral field of treatment in minimally invasive surgery and in aspirating particles. The plastic aspirator tips allow extremely accurate aspiration as the shape and design have been specially developed for working under a microscope or loupe. The tips are ideal because of their hygienic handling and ease of use.


- The green tip provides excellent contrast for precise aspiration
- Open field of view due to the narrow shape and conical tip
- Hygienic and time-saving, being sterile, individually packaged and ready for immediate use
- Various aspirator diameters for a wide range of treatments

## roeko Surgitip-endo

Surgitip-endo is an aspirating tip specially designed to dry root canals. The innovative design of the suction tip allows for a multi-directional flexibility, so it can easily be introduced into hard-to-reach root canals without having to bend the Canal Tip. Aspiration is highly effective at all angles and visibility is not constricted.

Surgitip-endo removes rinsing solutions and moisture quickly and effectively, shortening the drying time considerably. Less paper points are required. However, the use of paper points is recommended to check that the canal is thoroughly dry.

- Specially designed tip for multi-directional flexibility
- Effective aspiration and good canal accessibility
- Canal Tip ISO 60, inner diameter 0.35 mm
- Sterile, individually packaged and ready for use

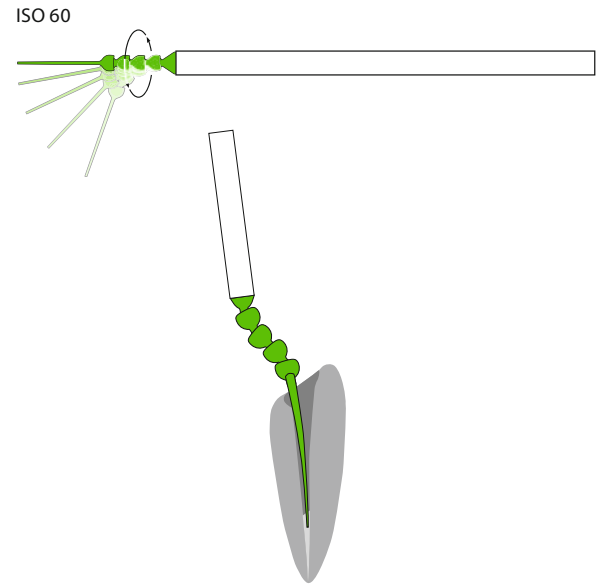


## roeko Surgitip-micro

With its 1.2 mm tip opening, the Surgitip-micro provides aspiration and precision during microsurgical procedures. The tip offers a good grip and is easy to manoeuvre, also making targeted aspiration possible in restricted dental operating environments. Surgitip-micro is sterile, individually packaged and ready for immediate use.




## roeko Surgitip

Surgitip is used for surgical procedures that require greater aspiration capabilities. It is available with 2.5 mm and 4.8 mm tip openings. The 4.8 mm tip consists of the bent aspirator tube without the green tip. Surgitip is sterile, individually packaged and ready for immediate use.

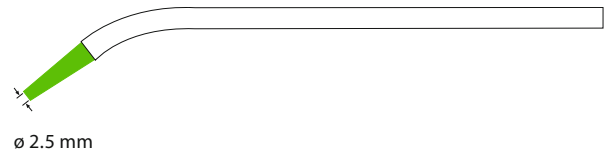


## roeko Aspirator tips

The non-sterile plastic aspirator tips are ideal for aspiration during non-surgical procedures. The length and shape correspond to those of the sterile 4.8 mm Surgitip. The aspirator tips are available in bulk packs with 100 pieces.

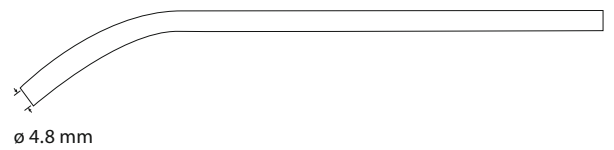


## roeko Double Adapter

The Double Adapter can be used to connect any ROEKO aspirator tip to standard aspirator systems with  $\varnothing$  11 mm or  $\varnothing$  16 mm tip openings. The Double Adapter is reusable and can be sterilised in the autoclave or chemiclave at up to 134 °C.

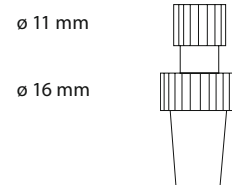

## roeko Concorde II

We recommend the Concorde II for universal aspiration. The plastic aspirator tip gently keeps the tongue and cheek away while aspirating large volumes. The tip made of high quality plastic is reusable and can be sterilised in the autoclave or chemiclave at up to 134 °C.

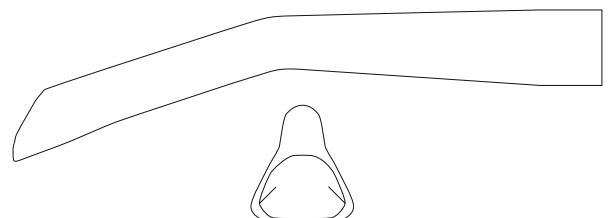


## ORDER INFORMATION

|                       | Tip opening ø mm | End opening ø mm | Contents               | REF                |
|-----------------------|------------------|------------------|------------------------|--------------------|
| <b>Surgitip-endo</b>  | 0.35 mm          | 6.5 mm           | 20 pcs *               | 6001 8280          |
| <b>Surgitip-micro</b> | 1.2 mm           | 6.5 mm           | 20 pcs *<br>100 pcs ** | 462 012<br>462 013 |
| <b>Surgitip</b>       | 2.5 mm           | 6.5 mm           | 20 pcs *<br>100 pcs ** | 462 025<br>462 026 |
|                       | 4.8 mm           | 6.5 mm           | 20 pcs *<br>100 pcs ** | 462 048<br>462 049 |
| <b>Aspirator tips</b> | 4.8 mm           | 6.5 mm           | 100 pcs *              | 461 048            |
| <b>Concorde</b>       |                  | 11 mm<br>16 mm   | 10 pcs ***<br>10 pcs   | 450 211<br>450 216 |
| <b>Adapter</b>        |                  | 11 mm<br>16 mm   | 1<br>1                 | 463 025<br>463 011 |
| <b>Double Adapter</b> |                  | 11 + 16 mm       | 1 Set                  | 463 000            |

\* incl. 1 Double Adapter | \*\*incl. 5 Double Adapters | \*\*\*incl. 1 Adapter



© COLTENE – www.coltene.com



**Coltène/Whaledent Ltd.**  
The President Suite  
Kendal House, Victoria Way  
Burgess Hill, West Sussex  
RH15 9NF / UK  
Tel +44 1444 235486  
Fax +44 1444 870640  
info.uk@coltene.com

**Coltène/Whaledent GmbH + Co. KG**  
Raiffeisenstraße 30  
89129 Langenau / Germany  
Tel +49 7345 805 0  
Fax +49 7345 805 201  
info.de@coltene.com

**Coltène/Whaledent Inc.**  
235 Ascot Parkway  
Cuyahoga Falls, Ohio 44223/USA  
Tel. + 1 330 916 8800  
Fax + 1 330 916 7077  
info.us@coltene.com

**COLTENE**