ultramat 2

HIGH SPEED AMALGAMATOR

HOCHFREQUENZ MISCHGERÄT AMALGAMADOR DE ALTA VELOCIDADE AMALGAMADOR DE ALTA VELOCIDAD **VIBREUR ELECTRONIQUE A GRANDE VITESSE** VIBRATORE PER AMALGAMA AD ALTA VELOCITÀ **HIGH SPEED CAPSULESCHUDDER HØJHASTIGHEDS-KAPSELRYSTER HØYHASTIGHETS-KAPSELRYSTER** NOPEA KAPSELISEKOITIN ΥΨΗΛΗΣ ΤΑΧΥΤΗΤΑΣ ΔΟΝΗΤΗΣ ΑΜΑΛΓΑΜΑΤΟΣ WYSOKOOBROTOWY MIESZALNIK NAGY SEBESSÉGŰ KAPSZULARÁZÓGÉP SUURE KIIRUSEGA AMALGAAMISEGISTI HITRI MEŠALEC AMALGAMA AUGSTAS FREKVENCES AMALGAMATORS DIDELIŲ APSUKŲ AMALGAMŲ MAIŠYTUVAS VYSOKORYCHLOSTNÍ AMALGAMÁTOR HITRI MEŠALEC AMALGAMA 高速アマルガムミキサー 高速银汞调拌机









5 year warranty amalgamator



- Safety feature
- Consistent accurate trituration
- Practical

L t a t a t

- Multi-use mixer
- High energy output
- Low energy usage
- 5 year warranty
- Simple to use
- Quiet low vibration
- Dual voltage

5 years warranty

Ultramat 2					
					5
Silamat S5*					
					3
AutoMix*					
					1
OptiMix*					
					1
ProMix*					
					1
RotoMix*					
					1
Dentomat Co	mpact*				
					1
Wig-L-Bug*					
					1
0 1	2	3	4	5	Years
	-	2	·	5	

ultramat 2 features:


Consistent and accurate trituration

The Ultramat 2 features a 0-16 second microprocessor-crystal controlled oscillator timer.



Multi-use mixer

The multi-use Ultramat 2 is suitable for all types of encapsulated dental materials.

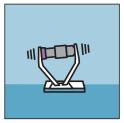

Quiet low vibration

The Ultramat 2 is the quietest, low vibration mixer.




Safety feature

The micro-switch automatically halts trituration when the lid of the Ultramat 2 is lifted.



Simple to use

The touch control system of the Ultramat 2 is quick and simple to operate.

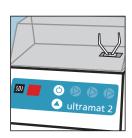

High energy output

The high energy Ultramat 2 has a mixing frequency of 4600 oscillations per minute +/-1%.



Practical

The one surface cover facilitates infection control cleaning and the non-porous plastic surface can be cleaned with any dental disinfectant.



Low energy usage The Ultramat 2 has economical low energy usage.



5 year warranty

The Ultramat 2's robust and reliable construction assures your purchase is a long lasting investment.



Dual voltage

The Ultramat 2 operates on voltages between 120 and 240v, 50 and 60Hz, for compatibility with multiple power voltages and frequencies.



t a g



Dimensions:

width	depth	height
217mm	173mm	180mm
8 1/2"	6 3/4"	7 1/8″

Voltage:

120 / 240 volts 50 / 60 Hz

Weight:

2.765kg 6.095lbs

Ultramat 2 (USA - includes power cord for USA) Reorder 5546045

Ultramat 2 (All countries - includes power cord for region product is shipped to) Reorder 5546052





DENTAL EQUIPMENT





*AutoMix, Dentomat Compact, OptiMix, ProMix, RotoMix, Silamat S5, Wig-L-Bug, are not registered trademarks of SDI Limited.

