







- 1 x VALO Grand LED curing light
- 4 x Batteries (rechargeable)
- 1 x Battery charger
- 1 x Power supply with universal plugs for battery charger 1 x Blue light blocking light shield
- 1 x Handpiece bracket holder
- 1 x Sample pack barrier sleeves
- **5972** Black
- **4865** Red Rock **4866** Midnight

• 4864 Sapphire

VALO Grand Corded Kit

- 1 x VALO LED curing light 7' cord 1 x Power supply with universal plugs 6' cord
- 1 x Handpiece bracket holder
- 1 x Blue light blocking light shield
- 1 x Sample pack barrier sleeves
- **5971** Black

VALO Cordless Kit 1 x VALO Cordless LED curing light

- 4 x Batteries (rechargeable)
- 1 x Battery charger
- 1 x Power supply with universal plugs for battery charger 1 x Blue light blocking light shield
- 1 x Handpiece bracket holder
- 1 x Sample pack barrier sleeves
- **5941** Black

VALO Corded Kit

- 1 x VALO LED curing light 7' cord 1 x Power supply with universal plugs 6' cord
- 1 x Handpiece bracket holder
- 1 x Blue light blocking light shield
- 1 x Sample pack barrier sleeves
- **5919** Black

VALO Grand Lens Accessories

4082 VALO Grand PointCure 2pk

5934 VALO PointCure 2pk

4081 VALO Grand ProxiCure 2pk

5936 VALO ProxiCure 2pk

4084 VALO Grand Translume *2pk*

5937 VALO Translume *2pk*

4319 VALO Grand Black Light 1pk

5939 VALO Black Light 1pk

4658 VALO Grand Interproximal *2pk*

4629 VALO Interproximal *2pk*

4628 VALO Grand White Light 2pk





ALL THE BENEFITS OF

THE REVOLUTIONARY, AWARD-WINNING

VALO™ CURING LIGHT,

WITH A 50% BIGGER LENS.

AVAILABLE ON ALL VALO CURING LIGHTS



UNPARALLELED POWER

Ultra-high-energy broadband LEDs produce a high-intensity light with an irradiance of 800–2300 mW/cm²



MULTIPLE WAVELENGTHS

Custom LEDs in three wavelengths polymerize all light-cured dental materials



DURABLE UNIBODY CONSTRUCTION

A bar of tempered, high-grade aerospace aluminum is CNC precision milled to create an extremely durable, lightweight handpiece



ACCESSIBILITY

The 11 mm head height allows easy and direct access to all curing sites



COLLIMATED BEAM

Optimally collimated beam delivers consistent, uniform power



MORE UNIFORM CURE

Unique glass lens system forms the collimated blended beam to give an even cure



AWARD-WINNING

The VALO family has won multiple Reality 5 Star Awards, Townie Choice Awards, and Pride Institute "Best of Class" Technology Awards, and has been named one of Dentistry Today's Top 50 Technology Products multiple times

EXCLUSIVE TO THE VALO GRAND LIGHT



LARGER TIP DIAMETER

The VALO Grand light has a 12 mm diameter lens, capable of covering 107 mm²



ADDITIONAL CURING BUTTON

Curing button on the underside facilitates intuitive operation



All VALO LED curing lights use a custom, multiwavelength light-emitting diode (LED) for producing high-intensity light at 385–515 nm, which is capable of polymerizing all light-cured dental materials.

Every member of the VALO curing light family also offers three curing modes to accommodate your preferences.

VALO GRAND LENS ACCESSORIES ADAPTABLE FUNCTIONS IN A SNAP





INTERPROXIMAL LENS

A focused 1 mm white light to help view interproximal caries.



TRANSLUME™ LENS

Aids in detection of cracks and defects in tooth or prostheses.



PROXICURE™ BALL LENS

Holds matrix against contact while curing composite for Class II restorations.



BLACK LIGHT LENS

Fluoresces chemicals in various dental materials.



POINTCURE™ LENS

Creates a 2.5 mm diameter curing area to aid in positioning and tack curing.





WHITE LIGHT LENS


Provides a pure natural light for accurate shade matching

AREA COVERAGE COMPARISON

10 MM MOLAR COVERAGE

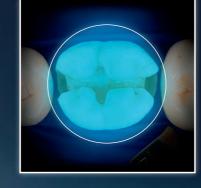

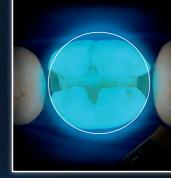

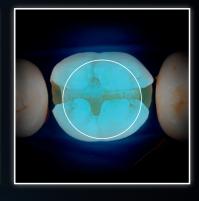





COMPETITOR

SURFACE AREA 107 mm² **78** mm²

46 mm²