


Cold-Curing acrylics

Vertex™ Castavaria, Vertex™ Castaquick & Vertex™ Castapress



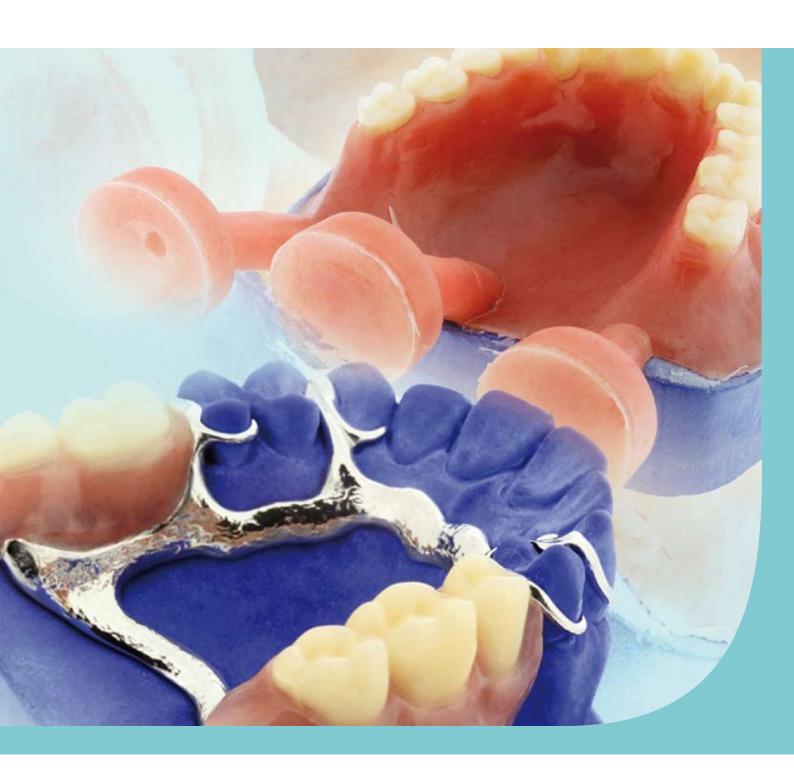

Cold-Curing acrylics

Cold-Curing for full & partial denture acrylics indications:

- Full and partial dentures
- Finishing of framework
- Repair, relining and rebasing of dentures
- Extensions to dentures

Vertex™ Castavaria



Vertex™ Castavaria is a multifunctional self polymerizing acrylic which is perfectly useable as a pouring, relining, rebasing and as a repair acrylic. The advantages of this acrylic are: minimized shrinkage, colour stable in the long run, stable polymerisation cycle with a perfect end result and the acrylic is long pourable. Besides Vertex™ Castavaria is modelable for a long period of time.

Vertex™ Castaquick



Vertex™ Castaquick provides more rapid settings, forms dough more quickly and is very suitable for adding artificial elements. Castaquick is only powder which must be used with the Castavaria liquid. For rebasing / re-lining dentures which are made by Castavaria and Castapress, Castaquick must be used.

Vertex™ Castapress



Vertex™ Castapress is a self polymerizing pouring/casting type denture base material also suitable for repair, relining, rebasing and extensions of partial dentures. The colour stability of the material is excellent because of the use of an unique accelerator system. Also available in Castapress Crystal Clear.

The following products are related to the Cold-Curing acrylics:



Vertex™ Casting Flask Using in combination with Vertex™ Castagel clear it is easy to control the whole process due to the transparency of the flask and the gel.



Vertex™ Castasil 21
Is a pour type of silicone meant for duplication of models.



Vertex™ Castagel
Can be used in
combination with Vertex™
Castavaria or Vertex™
Castapress to produce
effective a full or partial
denture.



Vertex™ Pouring cups
With the mixing cups for
powder and liquid, it is
possible to work with the
correct mixing ratio's easier.
Available in 100ml, 250ml
and 500ml.

Vertex-Dental B.V.

Johan van Oldenbarneveltlaan 62 | 3705 HJ Zeist P.O. box 10 | 3700 AA Zeist | The Netherlands Tel. +31 (0)30-69 76 749 | Fax +31 (0)30-69 55 188 info@vertex-dental.com | www.vertex-dental.com