**VITAPAN®** LINGOFORM

For all concepts of occlusion, thanks to multifunctional occlusal surface design.



# > FACTS & FEATURES

What?

• Natural, multifunctional premium posterior teeth made of MRP composite

### Benefits?

Setup is easy and intuitive

with "automatic" centric finding, thanks to occlusal surfaces designed according to the cogwheel principle

- Can be used universally and reliably for all concepts of occlusion thanks to multifunctional occlusal surface design
- Time-saving setup without extensive grinding, thanks to precisely defined centric
- Excellent durability thanks to high load capacity, abrasion stability and centric-supported function
- Good processability due to optimum material homogeneity and edge stability

What for?

• ideal solution for highly esthetic full, partial and implant-retained dentures.

# **VITAPAN® LINGOFORM**



### > BASIC INFORMATION

### **Characteristics**

#### Mould

- three-dimensional tooth geometry with natural lingobuccal dimension
- multifunctional occlusal surface design with defined centric contacts

### **Special features**

- finely nuanced surface texture for natural play of light
- Maxillary cusp angle: approx. 15° to 20°

### **Recommended indications?**

### Ideal prosthetic solution for:

- full and partial dentures
- implant-retained dentures

### Occlusal/setup concept

Suitable for all established concepts of lingual and balanced occlusion:

e.g., according to Gysi, Prof. Dr. A. Gerber, APF/TiF and ABC



lingualized occlusal concepts

### **Available shades**

- 27 x VITA SYSTEM 3D-MASTER® shades 0M3.1M1-5M3
- 15 x VITA classical A1–D4® shades A1–D4 (without B1)

### Available moulds/sets

- 4 x upper posteriors tooth moulds
- 4 x lower posterior tooth moulds

# > **OVERVIEW OF SHAPES**



\*) Bleached shade 0M3 available exclusively for selected mould(s).

# **VITAPAN® LINGOFORM**

# > RECOMMENDED COMBINATIONS

## **Combinations according to indications:**

## For full/partial dentures and implant-retained dentures:

- VITAPAN LINGOFORM Posterior with VITAPAN EXCELL Anterior
- See combination table for VITAPAN EXCELL Anterior

#### or

- VITAPAN LINGOFORM Posterior with VITAPAN PLUS Anterior
- See combination table for VITAPAN PLUS Anterior



## > ASSORTMENTS AND STORAGE





### **VITAPAN** LINGOFORM mould chart

Handy metal case with all moulds for tooth selection:

VITAPAN LINGOFORM Posterior


## VITAPAN LINGOFORM tooth storage

- X-small\*: 50 Sets / 400 Posteriors
- Small\*: 100 Sets / 800 Posteriors
- Medium\*: 200 Sets / 1,600 Posteriors
- Large\*: 300 Sets / 2,400 Posteriors

Available in VITA classical A1–D4® and VITA SYSTEM 3D-MASTER® shades \*\*

- \*) Note: You will receive the most frequently used moulds and shades with the packaged tooth storage. We are also happy to put together storage systems according to your individual specifications.
- \*\*) see available shades

# > SYSTEM COMPONENTS

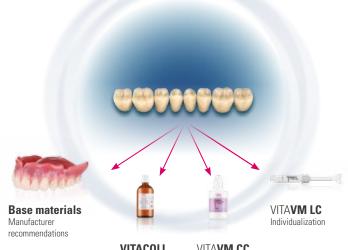





three-dimensional, natural moulds